# Optica<sup>™</sup> MDV & WDS

Illuminated Magnifier Examination Lights

# 

## **Product Description**

## Optica<sup>™</sup> MDV & WDS

Optica MDV Flexible LED Magnifier with Optimal Colour

**Optica WDS** Flexible Woods Light Magnifier with Optimal Colour for dermatology and veterinary use.



#### Product Images





Optica MDV

Optica WDS



## **Product Features**

#### Large Magnification

The large magnifier diameter of 157mm combined with 1.9x (3.5dpt) magnification allows the user to focus on fine details through a large-scale view.

# **Optimal Colour Rendering Properties**

6,500°K LEDs ensure less shadow formation and optimal colour rendering.

#### Precise Colour Display

The exceptional Ra of >93 accurately displays the high quality of white light and the R9 value >90 describes how accurately the user can differentiate between red tissues. CRI Ra >93. R9 >90. R13 >95

#### **Flexible Positioning**

Reinforced friction joints and a spring balanced arm allows a broad reach of up to 142cm without the lamp head drifting. This allows smooth handling for users to examine typically harder to reach places with ease.

#### Woods Light

The Woods light uses UV to pick up on skin conditions that may be missed using a standard LED light. Optica WDS can be used to diagnose a wide range of conditions such as;

- Favus
- Types of intertrigo
- Erythrasma

The powerful UV LEDs allow Optica WDS to be used in ambient light conditions without the need for darkened rooms.

#### Accurate Diagnosis

The Woods light uses UV to pick up on skin conditions that may be missed using a standard LED light.

#### **Visible Results**

Other Woods lights require the room to be completely darkened to use, yet due to powerful BLB LEDs within the WDS light, it can function without this.

#### **Approved Medical Device**

Optica is an approved medical device, which means that it can be used for a wide range of applications, including:

- Dermatology
- Human medicine
- Dentistry
- Veterinary medicine
- Monitoring sterilized products

# Optica<sup>™</sup> MDV & WDS

Illuminated Magnifier Examination Lights



Model Number MDV	Description		
OPTICAMDVT1C	Optica Medical Magnifier Ceiling Mounted Examination Lamp		
OPTICAMDVD	Optica Medical Magnifier Desk Mounted Examination Lamp		
OPTICAMDVW	Optica Medical Magnifier Wall Mounted Examination Lamp		
OPTICAMDVM	Optica Medical Magnifier Mobile Mounted Examination Lamp		

Model Number WDS	Description		
OPTICAWDST1C	Optica Woods Medical Magnifier Ceiling Mounted Examination Lamp		
OPTICAWDSD	Optica Woods Medical Magnifier Desk Mounted Examination Lamp		
OPTICAWDSM	Optica Woods Medical Magnifier Mobile Mounted Examination Lamp		
OPTICAWDSW	Optica Woods Medical Magnifier Wall Mounted Examination Lamp		

## **Technical & Performance Information:**

Unless otherwise stated: All linear dimensions in mm, angles in degrees (±5º) and weights in kg

## **OPTICA MDV**

General System Information	OPTICAMDVT1C	OPTICAMDVD	OPTICAMDVM	OPTICAMDVW		
Mount Type	Ceiling Mount	Desk Mount	Mobile Mount	Wall Mount		
Lamphead/Arm	OpticaMDVT1	OpticaMDV				
System Weight (kg)	8.2	2.5	9.6	2.2		
Conformity	CE (EU MDR 2017/745)					
Type of Luminaire	Examination (luminaire for diagnosis)					
Medical Device Classification	I					
Electrical Classification	I	I	ļ	I		
Electrical Conformity	BS EN 60601-1; BS EN 60601-2-41					
Additional Information	Permanent Electrical Connection Required	Power Lead supplied loose with UK plug				

#### **OPTICA WDS**

General System Information	OPTICAWDST1C	OPTICAWDSD	OPTICAWDSM	OPTICAWDSW		
Mount Type	Ceiling Mount	Desk Mount	Mobile Mount	Wall Mount		
Lamphead/Arm	OpticaWDST1	OpticaWDS				
System Weight (kg)	8.2	2.5 9.6		2.2		
Conformity	CE (93/42/EEC) including amendments 2007/47/EC					
Type of Luminaire	Examination (luminaire for diagnosis)					
Medical Device Classification	I					
Electrical Classification	I	l	ļ	I		
Electrical Conformity	BS EN 60601-1; BS EN 60601-2-41					
Additional Information	Permanent Electrical Connection Required	Power Lead supplied loose with UK plug				

Mount	Ceiling Mount	Desk Mount	Mobile Mount	Wall Mount
Mount Size HxWxD (mm)	40 to 166 x Ø310	150 to 210 x 70 x 40	1000 x Ø700	91 x 59 x 67
Weight of Mount (kg)	6.1	0.35	7.5	0.1
Fixing Type	Permanent	Moveable	Mobile	Permanent
Available Rotation Lateral about Spindle (±5°)	360	360	360	180
Interface Orientation	Vertical Down	Vertical Up	Vertical Up	Vertical Up
Colour	White	Black	Grey/White	RAL 7040- Telegrey 4
Vaterial	Steel	Polyamide / Steel	MIXED (Steel, plastic)	Plastic
Fastener Type	Ø16 mm Pin	Ø16 mm Pin	Ø16 mm Pin	Ø16 mm Pin
Additional information	Stem can be cut to length	-	-	-

See following page for Lamphead Technical Data

# Optica<sup>™</sup> MDV & WDS

Illuminated Magnifier Examination Lights


# 


Lamphead Photometric Data	OpticaMDV	OpticaWDS	
Description	LED Magnifying Examination Light with 'Daylight' white light	LED Magnifying Examination Light 'Daylight' white ligh and Woods light settings	
Central Illuminance @ 0.15 metre (Lux max ±10%)*	6,000	6,000	
Field Diameter D10 @ 0.15 meter (mm)	560	560	
Nominal Working Distance from Lamphead m	0.15	0.15	
Colour Rendition Ra (Av R1-R8)*	>93	>93	
Colour Rendition R9 (red)*	>90	>90	
Colour Rendition R13	98	98	
Colour Temperature (°K) *	6,500	6,500	
Number & Type of Lights per Lamphead	90 LED	90 LED	
A shatta a sh	Continuously dimmable	Continuously dimmable	
Additional	1.9x Magnification Acrylic Lens, 3.5 Dioptres	1.9x Magnification Acrylic Lens, 3.5 Dioptres	
LED Life (hours)	50,000	50,000	

\*-10% / +20%

Lamphead/Arm General & Mechanical Data	OpticaMDVT1	OpticaMDV	OpticaWDST1	OpticaWDS	
Arm Type	Articulated (Spring) for Ceiling Mounting	Articulated (Spring)	Articulated (Spring) for Ceiling Mounting	Articulated (Spring)	
Light On/Off Switch Location	Touch Panel on R	lear of Lamphead	Touch Panel on Rear of Lamphead		
Fastener Type	Ø16 m	ım Pin	Ø16 m	Ø16 mm Pin	
Base Joint Type and Angle (°)	Articulating 100°	Articulating 93°	Articulating 100°	Articulating 93°	
Base Rotation about Fastener (°)	360°	360° (max)	360°	360° (max)	
Elbow Angle (°)	177°	180°	177°	180°	
Arm Reach	1080		1080		
Lamphead Diameter (mm) Approx	220		220		
Lamphead Pitch (°)	180°		180°		
Lamphead Roll (°)	360°		360°		
Lamphead & Arm Weight (kg) Approx	2.1	2.1	2.1	2.1	
Power Supply Lead(s) [Plug]	Direct Supply	1.8m [UK]	Direct Supply	1.8m [UK]	

Lamphead/Arm Electrical Data	OpticaMDVT1	OpticaMDV	OpticaWDST1	OpticaWDS
Primary Power Supply	100-240VAC 50/60 Hz 100-240VAC 50/		AC 50/60 Hz	
Secondary Power Supply Unit	Integral 24V DC		Integral 24V DC	
Power Consumption @ Nominal Intensity	12W		12W	
Ingress Protection Rating	IP20		IP	20

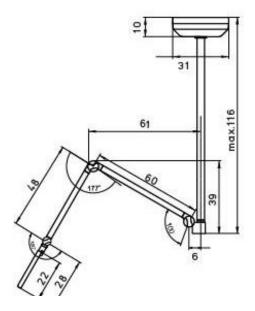

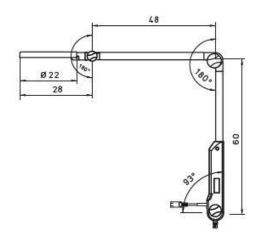

# www.curemednordic.se